


Human Apolipoprotein E III

APO-E3 recombinant protein	
Human Apolipoprotein E-III (Uniprot accession number P02649)	
Product Description	Recombinant human ApoEIII was expressed in <i>E. coli</i> with a purification tag at the N-terminus and purified using proprietary method including affinity, size exclusion chromatography. Purification tag was cleaved off by a specific protease.
Protein sequence	GKVEQAVETEPEPELRQQTEWQSGQRWELALGRFWDYLRWVQTLSEQVQE ELLSSQVTQELRALMDETMKELKAYKSELEEQLTPVAEETRARLSKELQAAQ ARLGADMEDVCGRLVQYRGEVQAMLGQSTEELRVRLASHLRKLRKLLRDA DDLQKRLAVYQAGAREGAERGLSAIRERLGPLVEQGRVRAATVGS LAGQPL QERAQAWGERLRARMEEMGSRTRDRLDEVKEQVAEVRACLEEQAAQIRLQ AEAFAQARLKSWEPLVEDMQRQWAGLVEKVQAAVGTSAAPVPSDNH
Predicted Molecular Mass	34.29 kDa
Expression system	<i>E. coli</i>
Purity	Greater than 90 % as determined by SDS-PAGE analysis
Labelling	No isotopic labelling
Tag information	Purification tag removed using a protease Glycine residue remains at the N-terminus
Protein Content	Quantitation is carried out by Bradford protein assay using BSA as standard.
Formulation	Lyophilized from 50 mM Tris-HCl, 300 mM NaCl, 1 mM DTT, pH 8.0 buffer. Centrifuge the vial prior to opening. Reconstitute by adding ultrapure water.
Shipping and storage	The product is supplied as a lyophilized powder and shipped at room temperature. Store at -80°C upon receipt. After reconstitution, protein can be kept at 4°C for a few weeks.
Usage	The product is for research use only. Not for diagnostic or therapeutic use.

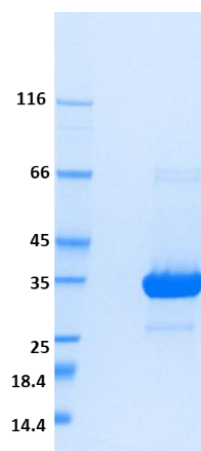

PROMISE Advanced Proteomics

Zone Minatec Entreprises – BHT 52A
7, parvis Louis Néel – CS20050
38040 Grenoble Cedex 9 – France

Tel : +33 (0) 4.38.02.36.50
Fax : +33 (0) 4.76.96.10.38
contact@promise-proteomics.com
www.promise-proteomics.com/