



Product name

DURVALUMAB Labelled standard

Catalog number IMU05501

Product description

Labelled Durvalumab is a recombinant, full length stable isotope labelled monoclonal antibody designed for use as an internal standard for quantitative analysis of Durvalumab in plasma samples. Durvalumab targets the immune checkpoint programmed deathligand 1 (PD-L1) and is used to treat a number of cancers

Similar products: Labelled Durvalumab-like, Labelled Nivolumablike, Labelled Ipilimumab-like, Labelled Pembrolizumab-like

mAb sequence

Heavy chain

EVQLVESGGGLVQPGGSLRLSCAASGFTFSRYWMSWVRQAPGKGLEWVANIKQDGSEKYYVDSVKGRFTISRDNAKNSLYLQMNSLRAE DTAVYYCAREGGWFGELAFDYWGQGTLVTVSSASTKGPSVFPLAPSSKSTSGGTAALGCLVKDYFPEPVTVSWNSGALTSGVHTFPAVLQS SGLYSLSSVVTVPSSSLGTQTYICNVNHKPSNTKVDKRVEPKSCDKTHTCPPCPAPEFEGGPSVFLFPPKPKDTLMISRTPEVTCVVVDVSHED PEVKFNWYVDGVEVHNAKTKPREEQYNSTYRVVSVLTVLHQDWLNGKEYKCKVSNKALPASIEKTISKAKGQPREPQVYTLPPSREEMTKN QVSLTCLVKGFYPSDIAVEWESNGQPENNYKTTPPVLDSDGSFFLYSKLTVDKSWQQGNVFSCSVMHEALHNHYTQKSLSLSPGK

Light Chain

EIVLTQSPGTLSLSPGERATLSCRASQRVSSSYLAWYQQKPGQAPRLLIYDASSRATGIPDRFSGSGSGTDFTLTISRLEPEDFAVYYCQQYGSL PWTFGQGTKVEIKRTVAAPSVFIFPPSDEQLKSGTASVVCLLNNFYPREAKVQWKVDNALQSGNSQESVTEQDSKDSTYSLSSTLTLSKADYE KHKVYACEVTHQGLSSPVTKSFNRGEC

Product features and protocols

Key features

Purity

>95%

as determined by SDS-PAGE

2

Labelling

 $Arg^{-13}C_{6}$, $^{15}N_4$ | Lys- $^{13}C_{6}$, $^{15}N_2$

Isotopic incorporation

as determined by LC-MS/MS

analysis of digested SIL-mAb

Other features

Expression System	Mammalian cells
Protein content	Quantitation is carried out by UV Absorbance at 280 nm
Formulation	5 mM Histidine pH 6 1% trehalose 0.002% Polysorbate 80

Template018-datasheet mab-V01

IMU05501_V01

Product preparation

For product preparation we recommend the following steps:


- 1. Briefly centrifuge the tube before opening
- 2. Reconstitute by adding the appropriate volume of ultrapure water for a final concentration of 200 μg/mL (e.g. 50 μl for 10 μg or 250 μl for 50 μg)
- 3. Vortex gently to insure complete dissolution
- 4. Wait 15 minutes at room temperature before proceeding further
- 5. Vortex gently again and centrifuge briefly

Product storage

The product is lyophilized and shipped at ambient temperature. Store at -80 °C upon receipt.

After reconstitution, the protein can be preserved at 4°C for a few weeks.

Avoid multiple freeze-thaw cycles

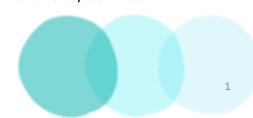





How to use our product



SIL mAbs are used as an internal standards by addition at the start of a sample preparation workflow. Both Labelled and target therapeutic mAbs will be processed throughout the entire analytical procedure controlling for variability caused by pre-analytical treatments or incomplete enzymatic digestion (1).

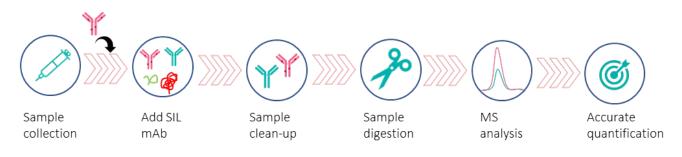

Supporting information

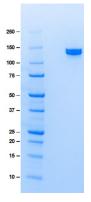

SDS-PAGE gel analysis of DURVALUMAB protein in non-reducing/ non heated conditions (NR/NH) Stained with Coomassie blue

References

1. D. Lebert, G.Picard, et al. Absolute and multiplex quantification of antibodies in serum using PSAQTM standards and LC-MS/MS, Bioanalysis. 2015, 7, 1237-51.



The product is intended for research use only. Not for diagnostic or therapeutic use.

Legal

This product is licensed under U.S. Patent No. 7,396,688 and foreign counterparts from E. I. du Pont de Nemours and Company. The purchase of this product conveys to the buyer the nontransferable right to use the purchased amount of the product for research and development only. Information about licenses for other uses is available from: E. I. du Pont de Nemours and Company; Attn: Associate Director, Commercial Development; DuPont Experimental Station E268; 200 Powdermill Rd.; Wilmington, DE 19803; 1-877-881-9787 (voice), 1-302-695-1437 (fax), licensing@dupont.com.