

Product description



Product name FREXALIMAB



Catalog number FRZ10201

FREXALIMAB is a second-generation monoclonal antibody that inhibits CD40 ligand.

It is developed for the treatment of relapsing multiple sclerosis.

mAb sequence

Heavy chain FREXALIMAB

EVQLQESGPGLVKPSETLSLTCTVSGDSITNGFWIWIRKPPGNKLEYMGYISYSGSTYYNPSLKRISISRDTSKNQFSLKSSVTAADTGVIYCAYSYGRTPYYFD
YWGQGTTLTVSSASTKGPSVFLAPSSKSTSGGTAALGCLVKDYFPEPTVSWNSGALTSVHFTPAVLQSSGLYSLSSVTVPSSSLGTQTYICNVNHKPSNTK
VDKKAEPKSCDKTHTCPPCPAPELGGPSVFLFPPKPKDTLMISRTPEVTCVVVDVSHRDPEVKFNWYVDGVEVHNAKTKPREEQYNSTYRVVSVLTVLHQD
WLNGKEYKCAVSNKALPAPIEKTISKAKGQPREPQVYTLPPSRDELTKNQVSLTCLVKGFYPSDIAVEWESNGQPENNYKTTTPVLDSGDSFFLYSKLTVDKSRW
QQGNVFCFSVMHEALHNHYTQKLSLSLSPGK

Light Chain FREXALIMAB

DIVMTQSPSFLSASVGDRTITCKASSNLGHAVAWYQQKPKGKSPKLLIYSASNRYTGVPDFRSGSGSGTDFTLTISLQPEDFADYFCQQYDDYPYTFGGGKLE
IKRTVAAPSVFIFPPSDEQLKSGTASVCLLNNFYPREAKVQWKVDNALQSGNSQESVTEQDSKSTYSLSSTLTLSKADYEKHKVYACEVTHQGLSSPVTKSFN
RGEN

Product features and protocols

Key features

- 1** Purity >95% as determined by SDS-PAGE
- 2** Labelling **Arg**-¹³C₆, ¹⁵N₄ | **Lys**-¹³C₆, ¹⁵N₂
- 3** Isotopic incorporation >98% as determined by LC-MS/MS analysis of digested SIL-mAb

Other features

Expression System	CHO cell line
Protein content	Quantitation is carried out by UV Absorbance at 280 nm
Formulation	5 mM Histidine pH 6 60 mM trehalose 0.01% Polysorbate 20

Product preparation

For product preparation we recommend the following steps:

- Briefly centrifuge the tube before opening
- Reconstitute by adding the appropriate volume of ultrapure water for a final concentration of 200 µg/ml (e.g. 50 µl for 10 µg or 250 µl for 50 µg)
- Vortex gently to insure complete dissolution
- Wait 15 minutes at Room temperature before proceeding further
- Vortex gently again and centrifuge briefly

Product storage

The product is lyophilized and shipped at room temperature.

Store at -80 °C upon receipt.

After reconstitution, the protein can be preserved at 4°C for a few weeks.

Avoid multiple freeze-thaw cycles


How to use our product



SIL mAbs allow to overcome the process variability since they are added at the very beginning of a sample preparation. This has positive impact on your analyte quantification, especially if the analyte interacts with other species commonly present within the matrix (1).



Supporting information

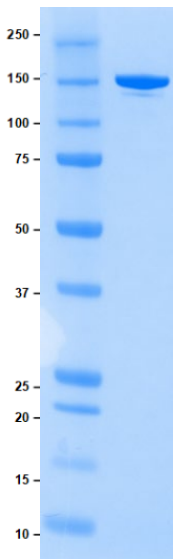

FREXALIMAB FRZ10201:
SDS-PAGE gel in
Reducing/Heated conditions
(RH) and non-reducing / non
heated conditions (NR/NH).
Stained with Coomassie blue.

References

1. **G.Picard, D. Lebert, et al.** PSAQ standards for accurate MS-based quantification of proteins: from the concept to biomedical applications, *J. Mass Spectrom.* 2012, 47, 1353-1363



The product is intended for research use only. Not for diagnostic or therapeutic use.

Legal

This product is licensed under U.S. Patent No. 7,396,688 and foreign counterparts from E. I. du Pont de Nemours and Company. The purchase of this product conveys to the buyer the nontransferable right to use the purchased amount of the product for research and development only. Information about licenses for other uses is available from: E. I. du Pont de Nemours and Company; Attn: Associate Director, Commercial Development; DuPont Experimental Station E268; 200 Powdermill Rd.; Wilmington, DE 19803; 1-877-881-9787 (voice), 1-302-695-1437 (fax), licensing@dupont.com.